
► Time-resolved Fluorescence Immunoassay (TRFIA) Method:

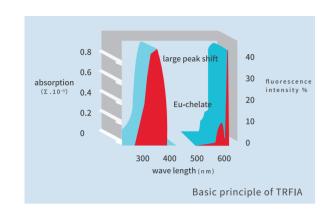

TRFIA is super-sensitive detection technique characterized by specific fluorescence of rare earth ions. It is not only highly sensitive, but also overcomes the instability of enzyme marker and is the better choice for immunological detection. The high fluorescence intensity and long life of labeled ionic chelates are beneficial to eliminate the influence of fluorescent substances in samples and environment on the test results.

▶ Analyzer Introduction:

The analyzer uses advanced methods of time-resolved fluorescence immunoassay (TRFIA) to perform in vitro quantitative detection of biomarkers of diabetes, infectious diseases, cardiovascular diseases, sex hormones, gastric diseases, thyroid function, tumors, COVID-19 and other diseases.

Application: laboratory, emergency, ICU, clinical laboratory, clinic, ambulance, field rescue and other clinical departments, etc.

► Application:



Ambulance








Test Items:

Category	Test Item
COVID-19	COVID-19 Antigen
	(COVID-19) IgM/IgG★
	SARS-CoV-2 Neutralizing Antibodies★
	COVID-19 Antigen (Self-test) (Nasal Swab)
	COVID-19 Antigen (Self-test) (Saliva)
	Inflfluenza A+B/COVID-19 Antigen
Diabetes	HbAlc★
	CRP★
	PCT★
Inflammation	SAA★
innammation	SAA/CRP ★
	PCT/CRP★
	IL-6
	cTnI★
	hs-cTnI
	СК-МВ
	Муо
Cardiac	cTnI/CK-MB/Myo ★
	NT-proBNP
	D-Dimer ★
	H-FABP
	BNP
	FDP★
	ST2
Daniel from aking	NGAL
Renal fun ction	mAlb

Category	Test Item
Thyroid	TSH★
	TT3
	TT4
	FT3
	FT4
	β-НСС+
	FSH
	PRL
Hormone	Progesterone
	LH
	АМН
Cashila Function	PG I /PG II
Gastric Function	G-17
Bone Metabolism	25-OH-VD
Anemia	Ferritin★
	B12
	\$100-β
Angiocarpy	НСҮ
	CCP★
Immuna System	IgG4
Immune System	RF★
	ASO★
Monkeypox Virus	Monleypox Virus Antibody(Latex)
	Monleypox Virus Antibody(Immunoassay)
	Monleypox Virus Antigen(Latex)

★Fingertip blood is available

New items are available soon!

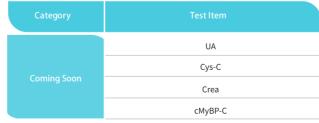

Category	Test Item
	сТпТ
	Testosterone
Coming Soon	E2

Lansion Biotechnology Co., Ltd.

Add: No.6 Qiande Road, Science Park, Jiangning District, 210000 Nanjing, Jiangsu Province, PEOPLE'S REPUBLIC OF CHINA

E-mail: biz@lansionbio.com
Web: en.lansionbio.com

Fax: +86-25-5875 8600

Tel: +86-25-5857 7600









PRODUCT CATALOGUE

POCT Point of Care Testing (Dry Fluorescence Immunoassay)



Quantitative, Sensitive, Reliable



Dry Fluorescence Immunoassay Analyzer (Portable)



▶ Features:

Method	Time-resolved Fluorescence Immunoassay (TRFIA)
Specimen	Serum/Plasma/Whole Blood/Nasal swab
Weight	1.3kg
Dimensions	225mm×152mm×105mm (L×W×H)
Screen	7 inch touch screen
Data Storage	10000
Printer	Built-in thermal printer
Battery	Built-in lithium battery
Communication	RS232(LIS/HIS), USB

Easy Operation:










Sample Dispense

QR Code Calibration Information Input

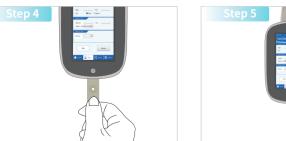

Automatic Detection

Result Output

LS-2100



Dry Fluorescence Immunoassay Analyzer (8 Channels)



- Sensitive
- Reliable
- High Throughput
- Rapid



▶ Features:

Method	Time-resolved Fluorescence Immunoassay (TRFIA)
Specimen	Serum/Plasma/Whole Blood/Nasal swab
Weight	7.5kg
Dimensions	390mm×290mm×205mm (L×W×H)
Screen	10.1 inch HD touch screen
System	Android
Test Speed	120T/H
Data Storage	10000
Printer	Built-in thermal printer
Communication	RS232(LIS/HIS), RJ45, USB, WIFI, Bluetooth

Easy Operation:





Insert Strip





Sample Dispense



Result Output

LS-4000



Dry Fluorescence Immunoassay Analyzer (Handheld)







Reliable

Sensitive

▶ Features:

Method	Time-resolved Fluorescence Immunoassay (TRFIA)
Specimen	Serum/Plasma/Whole Blood/Nasal swab
Weight	600g
Dimensions	191mm×84.5mm×45mm (L×W×H)
Screen	5.0 inch touch screen
Data Storage	50000
Battery	Built-in lithium battery (super standby time)
Communication	Type-C, 4G, USB, WIFI

(Up to 8 samples being detected simultaneously)